


# Titan Series Rescue Litters

11-0101A, Regular

11-0103A, Regular Split-Apart

11-0100A, Tapered

11-0104A, Tapered Split-Apart

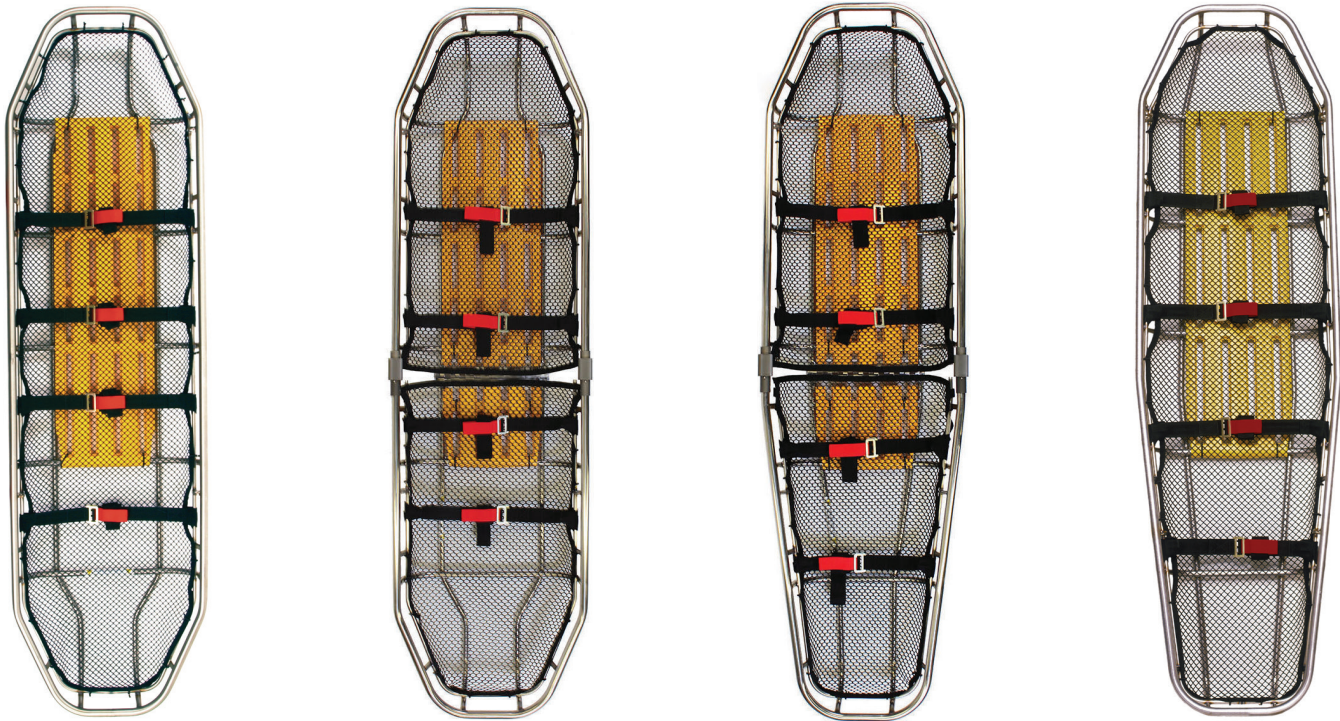

The **Titan Rescue Litter Series** is a complete rethinking of the classic Stokes military-style stretcher. The top rail is 1" (2.5 cm) diameter, making the stretcher easier to handle when carried for extended periods. Available in one-piece or two-piece and regular, or tapered models.

## Features

- Available in stainless steel
- Regular models are compatible with any backboard
- Patented StratLoad™ attachment points for secure attachment of lifting bridles
- Replaceable HDPE netting allows water and air flow plus a durable UHMW-PE back support
- Split-apart feature with LocSafe™ connectors provide fail-safe connection and allows the stretcher to be assembled in about 15 seconds
- Designed to surpass the US Military specification MIL-L-37957A. Tested to and beyond a 1,134 kg (2500 lb) static load test threshold
- Large 1" top tube for easy handling
- Titan Tapered enables easier access in confined areas
- Class 1 medical device

## Specifications

	Length	Width	Height	Weight	Load Rating
Regular	212.1 cm	59.4 cm	18.7 cm	14 kg	11 kn
	83.5 in	23.4 in	7.4 in	32 lb	2500 lb
Regular, Split-Apart	212.1 cm	61.6 cm	19.1 cm	15 kg	11 kn
	83.5 in	24.3 in	7.5 in	33 lb	2500 lb
Tapered	212.1 cm	59.4 cm	18.7 cm	14 kg	11 kn
	83.5 in	23.4 in	7.4 in	31 lb	2500 lb
Tapered Split-Apart	212.1 cm	61.6 cm	19.1 cm	15 kg	11 kn
	83.5 in	24.3 in	7.5 in	33 lb	2500 lb

