


Titan Titanium Rescue Litters

11-0116SL, Regular

11-0115SL, Regular Split-Apart

11-0118SL, Tapered

11-0119SL, Tapered Split-Apart

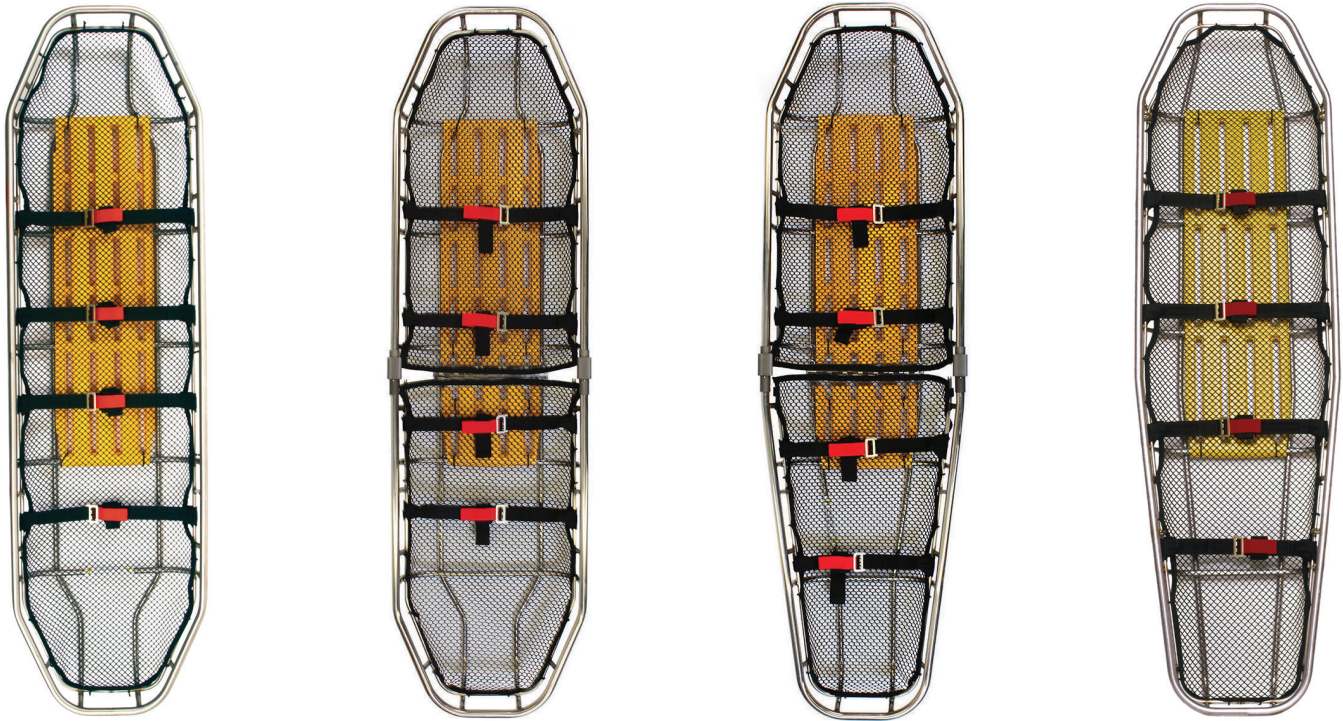

The **Titan-Ti Rescue Litter Series** is a complete rethinking of the classic Stokes military-style stretcher. The top rail is 1" (2.5 cm) diameter, making the stretcher easier to handle when carried for extended periods. Available in one-piece or split-apart and regular or tapered. For the ultimate in strength to weight ratio, Titan Titanium Rescue Litters set a new standard for rescue litters.

Features

- Titanium frame, lighter for bariatric use
- Regular model is compatible with any backboard
- Patented StratLoad™ attachment points for secure attachment of lifting bridles
- Replaceable HDPE netting allows water and air flow plus a durable UHMW-PE back support
- Split-apart feature with LocSafe™ connectors provide fail-safe connection and allows the stretcher to be assembled in about 15 seconds
- Designed to surpass the US Military specification MIL-L-37957A. Tested to and beyond a 1,134 kg (2500 lb) static load test threshold
- Large 1" top tube for easy handling
- Titan Tapered enables easier access in confined areas
- Class 1 medical device

Specifications

	Length	Width	Height	Weight	Load Rating
Regular	212.1 cm	59.1 cm	18.7 cm	8 kg	11 kn
	83.5 in	23.3 in	7.4 in	18 lb	2500 lb
Regular, Split-Apart	212.1 cm	60.3 cm	19.4 cm	9 kg	11 kn
	83.5 in	23.8 in	7.6 in	20 lb	2500 lb
Tapered	212.1 cm	59.1 cm	18.7 cm	8 kg	11 kn
	83.5 in	23.3 in	7.4 in	18 lb	2500 lb
Tapered, Split-Apart	212.1 cm	60.3 cm	19.4 cm	9 kg	11 kn
	83.5 in	23.8 in	7.6 in	19 lb	2500 lb

